




MANUAL PUNCH for Elmendorf and Bending

Code N3001

For easy cutting of samples for bending stiffness and tearing tests

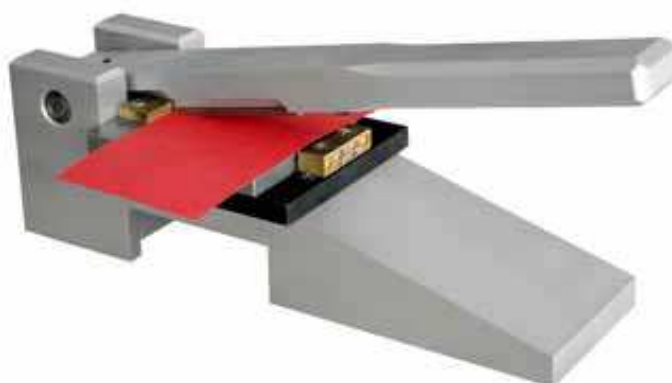

The material to be punched is placed on the lower blade

PTA SPAIN

C/ Gurutzegi N° 16 Pol. Belartza
20018 SAN SEBASTIÁN
(Guipúzcoa) - Spain

PTA FRANCE

14, rue Evariste Galois / ZAC les Ruieres
38320 EYBENS - France

PTA ASIA

50 East Coast Road #02-73 Roxy Square
428769 SINGAPORE - Singapore



DEVICE DESCRIPTION :

The hand punch allows fast preparation of samples for Elmendorf or bending tests. It consists of a triangular upper and lower blade. The punch is operated via a hand lever.

Choose dimensions when you order :

I x L : 38 x 75 mm
I x L : 38 x 70 mm
I x L : 50 x 62 mm
I x L : 50 x 65 mm
Other size on request

PROCESS DESCRIPTION :

The material to be punched is placed against the rest on the lower blade. When the hand lever is pushed down, the sample is evenly cut and falls through the die below the punch, where it can then be removed from the side.

WEIGHT AND DIMENSIONS (ONLY EQUIPMENT) :
100 x 300 x 200 mm (WxLxH) / 55 Kg

PTA SPAIN

C/ Gurutzegi N° 16 Pol. Belartza
20018 SAN SEBASTIÁN
(Guipúzcoa) - Spain

PTA FRANCE

14, rue Evariste Galois / ZAC les Ruieres
38320 EYBENS - France

PTA ASIA

50 East Coast Road #02-73 Roxy Square
428769 SINGAPORE - Singapore